


Maximized in Performance,
Minimized in Size

SE-301

Electrocardiograph

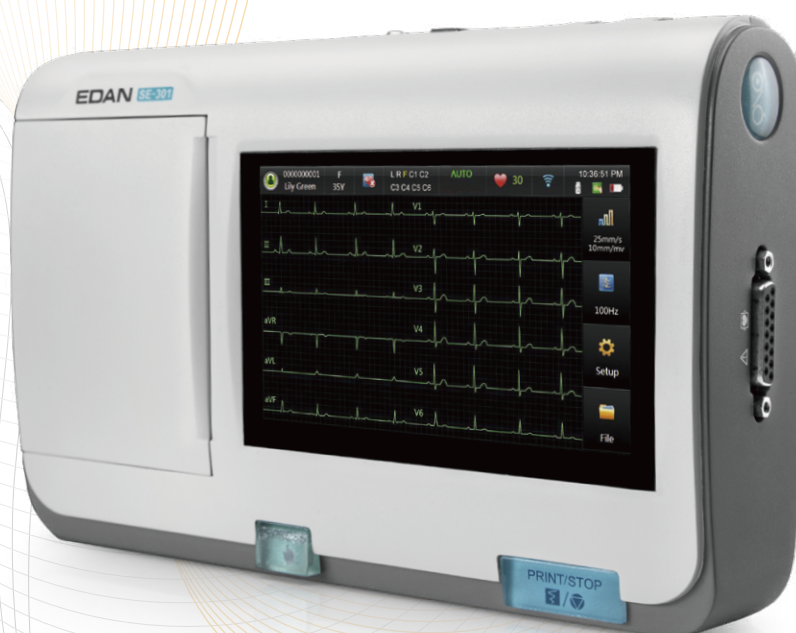

SE-301

3-Channel Electrocardiograph

SE-301 is the new generation 3-channel ECG, which incorporates the benefits of the advanced technology and innovative design in an extremely compact and lightweight device. It is specifically developed for those who pursue superb portability without compromising on the quality for mobile environments. A variety of fantastic features are embodied in an amazing small body, making it more than just a 3-Channel ECG.




Less than 1 kg



Full-touch
Operation



One-button
Operation

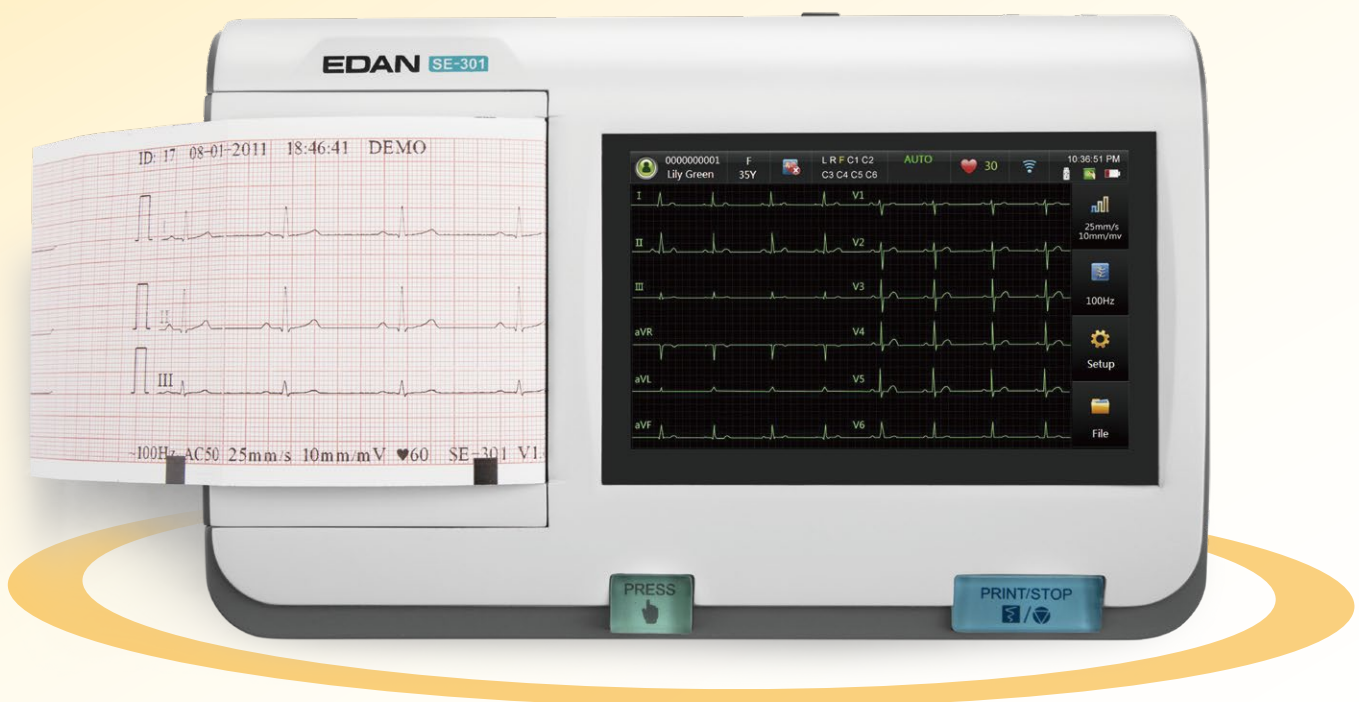

All for a reliable result

CMRR

A high Common Mode Rejection Ratio up to 140dB restrains signal noises to the greatest extent, thus offers you traces with excellent quality.

Bandwidth

A wide bandwidth from 0.01 to 300 Hz guarantees the accuracy of the ECG signal as well as the interpretation especially for pediatric patients.

Sampling Rate

The leading sampling rate as high as 16,000 Hz aims at capturing all the details to generate precise result.

Convenient Information Communication

Powerful External Device Support

SE-301 arms with two USB ports and one microSD card slot, which extends its capability by supporting various external devices such as barcode scanner, external printer, USB Memory Drive as well as microSD card.

Bi-directional Communication with Data Management System

Via the LAN, build-in Wi-Fi or built-in GSM/GPRS cellular modem, SE-301 is able to receive orders from SE-1515 Data Management System and transmit ECG reports back to have further analysis, editing, saving or printing.

Multi-format Report

SE-301 is able to export ECG reports in standard PDF format and optional SCP, XML and DICOM format, which offers you greater capability in integrating the ECGs into different data management systems.




Signal Quality Detection

By showing different colors to indicate normal, poor or missing signals, this function helps you to identify the signal quality even at a glance.

Reliable ECG Algorithm

SEMIP® ECG algorithm offers you an accurate and reliable reference to maximize clinicians' confidence in making diagnosis. Glasgow Algorithm with specialized Children ECG Interpretation is also optional for SE-301.

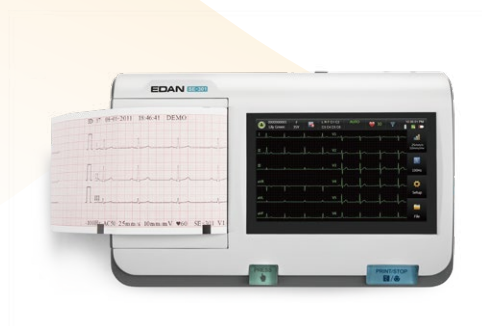

SE-301

Electrocardiograph

Global Headquarters:

Edan Instruments, Inc. | 15 Jinhui Road, Pingshan District, Shenzhen
518122 P.R. China | +86.755.26898326 | www.edan.com | info@edan.com

U.S. and Canada inquiries:

EDAN Diagnostics, Inc. | 9918 Via Pasar, San Diego, CA 92126
+1.858.750.3066 | www.edandiagnostics.com | edan-info@edandiagnostics.com

© Edan Instruments, Inc. All rights reserved. Features and specifications are subject to change without prior notice. No reproduction, copy or transmission may be made without written permission.
Not all products or features are available in all countries, contact Edan for local availability.